

Universität Zürich
Zentrum für Zahnmedizin
Klinik für Präventivzahnmedizin, Parodontologie und Kariologie (PPK)
Direktor: Prof. Dr. med. dent. Thomas Attin
Computergestützte Restaurative Zahnmedizin
Leitung: Prof. Dr. med. dent. Dr. rer. biol. hum. A. Mehl

Arbeit unter Leitung von Prof. Dr. med. dent. Dr. rer. biol. hum. A. Mehl

A 3-dimensional accuracy analysis of chairside CAD/CAM milling processes

INAUGURAL-DISSERTATION
zur Erlangung der Doktorwürde der Zahnmedizin
der Medizinischen Fakultät
der Universität Zürich

vorgelegt von
Gabriel Tobias Bosch
aus Deutschland

Genehmigt auf Antrag von Prof. Dr. med. dent. Dr. rer. biol. hum. A. Mehl
Zürich 2014

Publikationshinweis

A 3-dimensional accuracy analysis of chairside CAD/CAM milling processes

Akzeptiert am: 15.05.2014
Epub ahead of print: 30.06.2014

Journal: The Journal of Prosthetic Dentistry



A 3-DIMENSIONAL ACCURACY ANALYSIS OF CHAIRSIDE CAD/CAM MILLING PROCESSES

Gabriel Bosch, MedDent,^a Andreas Ender, DrMedDent,^b and Albert Mehl, Prof Dr, DrMedDent^c

Center of Dental Medicine, University of Zürich, Zürich, Switzerland

Statement of problem. Milling is a central and important aspect of computer-aided design and computer-aided manufacturing (CAD/CAM) technology. High milling accuracy reduces the time needed to adapt the workpiece and provides restorations with better longevity and esthetic appeal. The influence of different milling processes on the accuracy of milled restorations has not yet been reviewed.

Purpose. The purpose of this study was to investigate the influence of different milling processes on the accuracy of ceramic restorations.

Material and methods. Four groups of partial crowns were milled (each n=17): Three groups in a 4-axial milling unit: (1) 1-step mode and Step Bur 12S (12S), (2) 1-step mode and Step Bur 12 (1Step), (3) 2-step mode and Step Bur 12 (2Step), and (4) one group in a 5-axial milling unit (5axis). The milled occlusal and inner surfaces were scanned and superimposed over the digital data sets of calculated restorations with specialized difference analysis software. The trueness of each restoration and each group was measured. One-way ANOVA with a post hoc Tukey test was used to compare the data ($\alpha=.05$).

Results. The highest trueness for the inner surface was achieved in group 5axis (trueness, $41 \pm 15 \mu\text{m}$, $P<.05$). The 4-axial milling unit exhibited trueness at settings ranging from $61 \mu\text{m}$ (2Step) to $96 \mu\text{m}$ (12S). For the occlusal surface, the highest trueness was achieved with group 5axis (trueness, $42 \pm 10 \mu\text{m}$). The 4-axial milling unit exhibited trueness at settings ranging from $55 \mu\text{m}$ (1Step) to $76 \mu\text{m}$ (12S).

Conclusions. Restorations milled with a 5-axial milling unit have a higher trueness than those milled with a 4-axial milling unit. A rotary cutting instrument with a smaller diameter results in a more accurate milling process. The 2-step mode is not significantly better than the 1-step mode. (J Prosthet Dent 2014;■:■-■)

CLINICAL IMPLICATIONS

The goal of the milling process is to generate an exact copy of the digitally calculated restoration. The 5-axial milling unit came closer to the digitally calculated restoration than the 4-axial milling unit.

Dental ceramics have proven longevity and remain the materials of choice for esthetic restorations.¹⁻⁵ Ceramic restorations can be produced in different ways, one of which involves computer-aided design and computer-aided manufacturing (CAD/CAM). The goal of CAD/CAM technology is standardized, reproducible production that is both efficient and accurate.⁶

Clinical long-term success depends on factors such as adequate cementation, restoration design, preparation design, and (importantly) the marginal and internal fit of the restoration to the tooth. An important step with regard to the fit of the restoration is the fabrication process. Ceramic crowns can be distorted during this process, which can negatively affect the fit and compromise the success of the restoration.⁷⁻¹⁰

Important considerations are how accurate the milling process is and whether it damages the restoration.

CAD/CAM allows the dentist to work chairside and is the fastest way to produce individual ceramic restorations. Cerec (Sirona Dental Systems) is a well-known chairside CAD/CAM system that has proven both reliable and efficient.¹¹⁻¹³ The optical scanning system (Cerec Bluecam; Sirona Dental

^aResearch Assistant, Division for Computerized Restorative Dentistry.

^bResearch Assistant, Division for Computerized Restorative Dentistry.

^cHead of Division, Division for Computerized Restorative Dentistry.

Systems) has also been found to produce reliable and highly accurate 3-dimensional digital impressions.¹⁴⁻¹⁸ The adaption of the restoration to the preparation with regard to marginal and internal fit has been found to be clinically reliable in numerous studies.¹⁹⁻²³ However, marginal gap values range from 35 to 246 μm , and internal gaps range from 17 to 206 μm .²⁴ Relatively large internal gaps can result from the milling process,²⁵ but gaps of up to 150 μm are clinically tolerable.²⁶ The marginal gap differs among CAD/CAM systems,²⁷⁻²⁹ but with the continual development of milling machines, marginal gaps are becoming progressively smaller.³⁰

The accuracy of fitting to the underlying tooth structure is an essential consideration and could affect the longevity of restorations.^{31,32} Larger gaps are associated with accelerated plaque accumulation, secondary caries, marginal discoloration, exposure of the luting resin, dissolution of the cement, and increased risk of microleakage and microcracks.³³⁻³⁵ When the marginal gap is greater than 100 μm , removing excess cement is more difficult.³⁶ A restoration with inadequate fit may lead to marginal chipping, and even small chips can result in the later clinical failure of ceramic restorations.³⁷

With the CAD/CAM technique, high-precision scans are possible, and the software can calculate a restoration with exact control over contact points and design. The purposes of this study were (1) to evaluate the trueness of different milling processes that use 2 different milling units by comparing the milled restoration with the original digital data set and (2) to visualize the deviations and marginal chippings caused by the milling process. The primary null hypothesis was that no quantitative differences would be found between the 2 different milling devices investigated or their respective milling processes.

MATERIAL AND METHODS

Seventeen clinical ceramic preparations were selected for this study

(2 inlays with 2 surfaces, 5 inlays with 3 surfaces, 6 inlays with 4 surfaces, 4 onlays). All preparations were performed in accordance with the manufacturer's guidelines. The preparations were digitally copied by using an intraoral scanning system (Cerec Bluecam; Sirona) and were printed in acrylate polymer (Objet MED610; Objet Geometries GmbH) with a 3-dimensional printer (Objet Eden 260V; Objet Geometries GmbH).

The printed casts were lightly powdered (Cerec Optispray; Sirona Dental Systems), and all were scanned with a digital intraoral scanning system (Cerec Bluecam Connect, v4.03; Sirona Dental Systems) by the same experienced operator. The preparation border was defined, and then the data set was transferred to CAD/CAM software (inLab 4.0, v4.02; Sirona Dental Systems). The restorations were calculated and sent to the milling preview.

The 4 milling procedures investigated are displayed in Table I. Groups 12S, 1Step, and 2Step were milled with a 4-axial milling unit (inLab MC XL; Sirona Dental Systems). This milling unit uses 2 instruments for the milling process: a stepped bur (Step Bur), milling only the inner surface, and a pointed bur (Cylindrical Pointed Bur), milling only the outer surface of the restoration. The Step Bur is available in 2 diameters, Step Bur 12S (1.2 mm) and Step Bur 12 (1.0 mm). In group 12S, the restorations were milled with a Step Bur 12S and a Cylindrical Pointed Bur 12S. In group 1Step, the restorations were milled with a Step Bur 12 and a Cylindrical Pointed Bur 12S. In group 2Step, the restorations were milled with a Step Bur 12 and a Cylindrical Pointed Bur 12S by using a

so-called 2-step milling procedure. For those 3 groups, the parameters for the restorations were set as follows: Spacer, 80; Marginal-Adhesive Gap, 60; Occlusal-Milling Offset, 0; Proximal-Contacts Strength, -50; Occlusal-Contacts Strength, 25; Minimal Thickness (Radial), No; Minimal Thickness (Occlusal), No; Margin Thickness, 20; Consider Instrument Geometry, Yes; Remove Undercuts, Yes.

In group Arctica, the restorations were milled with a 5-axial milling unit (Arctica; KaVo). The restorations were exported to Surface Tessellation Language (STL) files (a standard format for CAD/CAM data exchange) with the manufacturer's recommended parameters as follows: Spacer, 20; Marginal-Adhesive Gap, 0; Occlusal-Milling Offset, 0; Proximal-Contacts Strength, -50; Occlusal-Contacts Strength, 25; Minimal Thickness (Radial), No; Minimal Thickness (Occlusal), No; Margin Thickness, 0; Consider Instrument Geometry, No; Remove Undercuts, Yes.

The STL files were then imported into the 5-axial CAM software and milled without any further changes. For all test groups, the calculated milling surfaces of the restorations were exported as STL files for later comparison with the milled restorations.

After milling, the sprue was removed from the restorations, and they were placed in the printed casts without any further adjustments. All milled restorations were scanned with a 3-dimensional scanning system (Bluecam, inLab 4.0, v4.02; Sirona Dental Systems). Bluecam is a scanner with a trueness for single-tooth scans of $\pm 19.2 \mu\text{m}$.¹⁷ The surfaces were lightly powdered (Cerec Optispray; Sirona Dental Systems), the restorations were scanned from the

TABLE I. Test groups with different settings (n=17)

Group	Rotary Instrument	Milling Unit	Milling Option
12S	Step Bur 12S	inLab MC XL	1-step
1Step	Step Bur 12	inLab MC XL	1-step
2Step	Step Bur 12	inLab MC XL	2-step
Arctica	Arctica Bur Set	Arctica	5-axis

occlusal view, and the 3-dimensional data sets were exported as STL files (occlusal surface comparison).

The milled restorations were then fixed to an object plate provided with reference grooves, with the inner surfaces of the restorations facing upward. The surfaces were lightly powdered (Cerec Optispray) and scanned, and the 3-dimensional data sets were exported as STL files (inner surface comparison).

Within each group, the scanned surfaces of the milled restorations were compared with the calculated milling surfaces with specialized difference analysis software (OraCheck v1.00.10; Cyfex). The software superimposes the STL files by using a best-fit algorithm for closest-point matching of both surfaces. The occlusal or inner surfaces were selected precisely and the restorations superimposed. The software calculated and measured the distance (positive or negative) from every surface point (approximately 20 000 per surface matching) from the milled to the original surface. The point-by-point difference values for each single superimposition were exported to an American Standard Code for Information Interchange (ASCII) text file (comma-separated values). All the calculated distances were imported into statistical software (IBM SPSS v19.0; IBM Corp). For each superimposition of 2 surfaces, the 10th and 90th percentiles were calculated. The metric value for the deviation between 2 surfaces was defined as the $(90\%-10\%)/2$ percentile (deviation measure [DM]). This DM gives the level by which approximately 80% of the matched area has less negative and positive deviations. In the next step, the mean value and the SD of the DMs were calculated for each group. These mean values describe the trueness of the milling process in terms of the deviation from the calculated original data. A lower value indicates a more accurate milling process. To compare the different groups, 1-way ANOVA with the post hoc Tukey honestly significant difference test was used as a statistical test ($\alpha=.05$).

All visual examinations of the differences between the original and the milled restorations were performed by the same experienced dentist by means of color-coded difference images with boundary values set from $+100 \mu\text{m}$ to $-100 \mu\text{m}$. Separate aspects of the images were examined, as shown in Table II. Aspect 1 was the occlusal relief, aspect 2 was the fine structures in the inner surface, aspect 3 was the surfaces with an angle close to the insertion axis in the inner surface, and aspect 4 was the marginal area in the inner surface up to $100 \mu\text{m}$ from the outer edge. The different aspects were visually evaluated separately for each restoration in direct comparisons between all groups. For the comparison of each restoration (4 images), scores were assigned as follows: 1 for the best

result, 4 for the worst result. Ranking was determined according to the amount of pink, green, and purple present. Pink indicates a loss of $\geq 100 \mu\text{m}$ of material, green indicates almost no deviation from the CAD data set, and purple indicates differences $\geq 100 \mu\text{m}$ larger than the original. The same score in different groups means no visual difference in quality. A ranking between the groups was calculated based on the sum of the restorations for each group.

RESULTS

The results of the visual examination are shown in Table III and the trueness results for all groups are displayed in Table IV. Figure 1 shows graphs of the means and SDs for trueness. Table V

TABLE II. Aspects for visual examination

Aspect	Description
Aspect 1	Occlusal relief
Aspect 2	Fine structures in inner surface
Aspect 3	Surfaces with angle close to insertion axis in inner surface
Aspect 4	Marginal areas in inner surface up to $100 \mu\text{m}$ from outer edge

TABLE III. Ranking of visual examinations


Ranking	Step	Step	Step	Arctica
	Bur 12S	Bur 12, 1-Step	Bur 12, 2-Step	
Aspect 1	4	2	2	1
Aspect 2	4	3	2	1
Aspect 3	4	3	2	1
Aspect 4	4	2	2	1

1 = best results; same number = no difference.

TABLE IV. Trueness for all test groups (trueness with $[90\% - 10\%]/2$ percentile)

Group	Process	Inner Surface	
		Occlusal Trueness	Trueness
12S	MC XL, Bur 12S, 1-step	$76 \mu\text{m} \pm 40 \mu\text{m}$	$96 \mu\text{m} \pm 68 \mu\text{m}$
1Step	MC XL, Bur 12, 1-step	$55 \mu\text{m} \pm 18 \mu\text{m}$	$67 \mu\text{m} \pm 20 \mu\text{m}$
2Step	MC XL, Bur 12, 2-step	$67 \mu\text{m} \pm 24 \mu\text{m}$	$61 \mu\text{m} \pm 22 \mu\text{m}$
Arctica	Arctica	$42 \mu\text{m} \pm 10 \mu\text{m}$	$41 \mu\text{m} \pm 15 \mu\text{m}$

MC XL, inLab MC XL (Sirona Dental Systems).

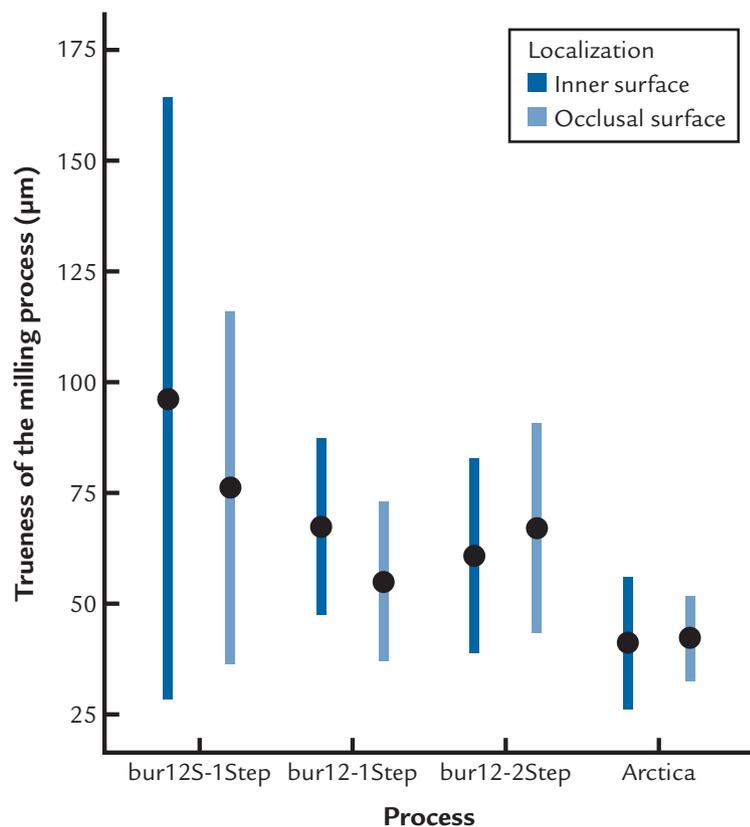

1 Trueness values and SDs for all test groups.

TABLE V. All test groups exhibiting statistically significant differences

Differences	P
Inner Surface	
Arctica versus Bur 12S, 1Step	<.001
Bur 12S versus Bur 12, 2Step	.040
Occlusal Surface	
Arctica versus Bur 12, 2Step	.029
Arctica versus Bur 12S, 1Step	.001

shows the groups that were statistically significantly different. For the occlusal area, df was 3 and F was 5.802. For the inner surface, df was 3 and F was 6.166. For all test groups, the trueness in the occlusal areas ranged from 42 to 76 μm ; in the inner areas, from 41 to 96 μm . Group Arctica exhibited the highest trueness in both areas, with 42 \pm 10 μm in the occlusal area and 41 \pm 15 μm in the inner surface. For the inner surface, group Arctica was statistically significantly better than group 12S ($P < .001$). For the occlusal surface, group Arctica was statistically significantly better than


groups 2Step ($P = .029$) and 12S ($P = .001$). The difference images from group Arctica showed the most accurate results in all 4 aspects (Figs. 2D, 3D). Group 2Step yielded a trueness of 67 \pm 24 μm for the occlusal area and 61 \pm 22 μm for the inner surface. In both areas, group 2Step did not differ statistically significantly from group 1Step, but in the inner surface, it was statistically significantly better than group 12S ($P = .040$). The difference images show the second-best results in all 4 aspects (see Figs. 2C, 3C). In aspect 1 (the occlusal relief) and in aspect 4 (the marginal areas), the results for groups 1Step and 2Step were the same. Group 1Step yielded a trueness of 55 \pm 18 μm in the occlusal area and 67 \pm 20 μm in the inner surface. The difference images were close to the difference images of group 2Step (see Figs. 2B, 3B). However, with regard to the fine structures in the inner surface in aspect 2 and to the surfaces with an angle close to the insertion axis in aspect 3, group 2Step was somewhat better. Group 12S exhibited the

lowest trueness in the occlusal area (76 \pm 40 μm) and in the inner surface (96 \pm 68 μm). The difference images shown in Figures 2A and 3A display the largest deviations compared with the CAD data set in all 4 aspects.

DISCUSSION

The trueness of chairside milling processes was examined in this in vitro study. On the basis of the results of this in vitro study, the null hypothesis, that no quantitative differences would be found between the 2 different milling devices investigated or their respective milling processes, has to be rejected. However, drawing meaningful comparisons with previously reported research is problematic. In previous studies, the marginal fit and internal fit were measured. Conventional methods measure these parameters in only 1 or 2 dimensions, and they measure the gap between the restoration and die with microscopy at 2 to 150 points in 5 to 10 specimens per group. Additionally, previous studies have evaluated fit with the elastomeric putty-wash technique, with low-viscosity silicone to duplicate the cement space and evaluate it photometrically or analyze it based on its density and weight.^{24,25}

Achieving clinically relevant results may require between 50 and 230 measuring points.³⁸ Modern computer-aided techniques can better evaluate the fit of the restoration, because they yield much more extensive and informative data in 3 dimensions.³⁹ Three-dimensional analysis has also proven valid and reliable.⁴⁰ In previous studies, the internal and marginal gaps have been examined with 3-dimensional analysis by superimposing scans of the die over scans of the inner surface from the ceramic crown, or by superimposing scans of the die over scans of the fit-checker on the die to measure the thickness of the fit-checker, which replicates the cement space.^{24,39} With an accurate intraoral scanning system, the fit of the restoration in the oral cavity depends on the milling process.¹⁸ Thus, this study did not focus directly on the marginal and internal fit

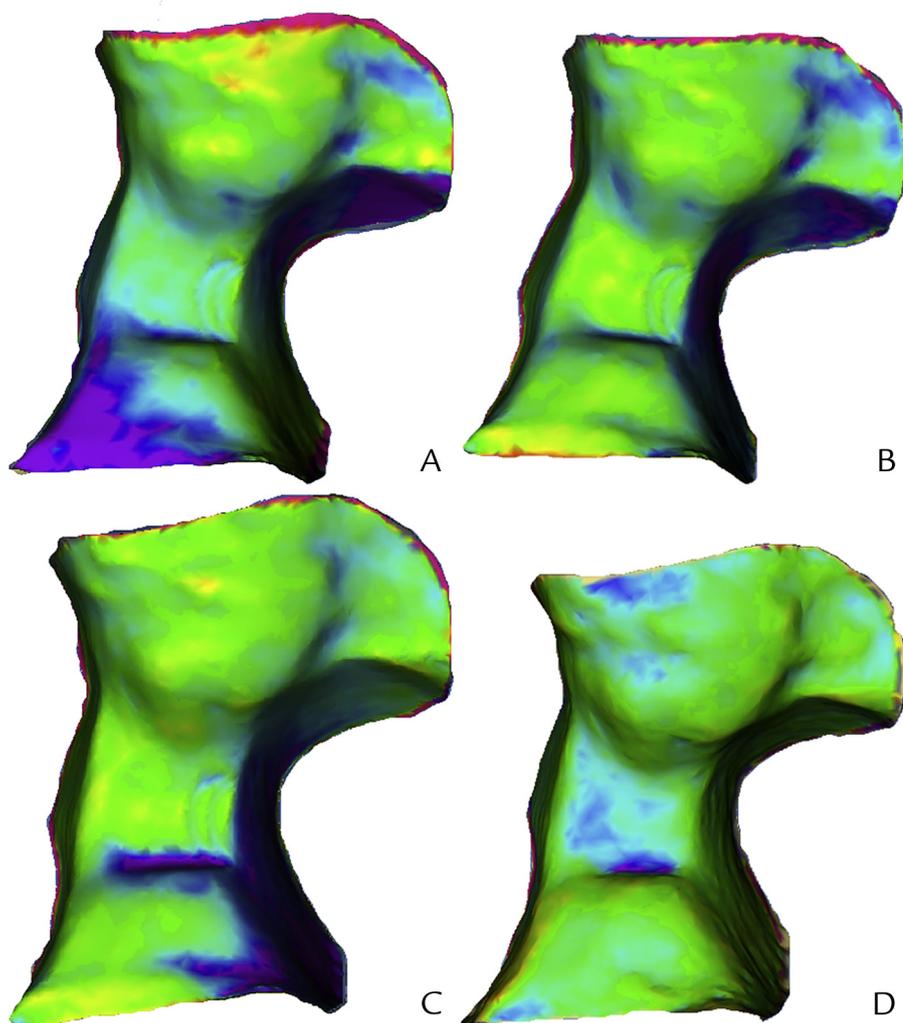

2 Deviations between digital calculated cast and scanned milled cast (trueness). A, Bur 12S, 1-step mode. B, Bur 12, 1-step mode. C, Bur 12, 2-step mode. D, Arctica. Color-coded from $-100\ \mu\text{m}$ (blue) to $+100\ \mu\text{m}$ (red).

but rather on the milling process. With better trueness in the milling process, a better fit is facilitated. An inaccurately milled restoration results in a poor fit, with numerous occlusal and proximal contacts, which evidently differs substantially from the contact points precisely determined by the CAD/CAM software. Correcting the contact points intraorally may negatively affect the esthetics of the restoration, the contact distribution, the chairside time, and the longevity.^{31,32}

In this study, the point-by-point differences were measured between the digitally acquired data and the milled restoration caused by the milling process. With the method reported, the inner and outer surfaces could be measured in addition to the


internal fit. Another advantage of this method is that unlike other methods, this one requires neither fit-checker nor several dies, making it easier and more fail-safe. The difference images facilitate direct visual feedback encompassing the entire restoration and make it possible to locate the more imprecise areas. Bluecam has been validated as being accurate for single-tooth scans, with a trueness of $\pm 19.2\ \mu\text{m}$ in this context, and is sufficiently accurate for use in this measuring method.¹⁷

The 4-axial milling unit (inLab MC XL; Sirona Dental Systems) uses the Step Bur for the inner area and the Cylinder Pointed Bur for the occlusal area. During milling, the rotation axis of the material holder is fixed to

ensure that both instruments maintain an appropriate insertion axis. In 2-Step mode, the restoration is milled in 2 cycles. In the first milling cycle, the restoration is milled with an additional material thickness of 0.3 mm at all surfaces. The second milling cycle removes the remaining 0.3 mm of material to get the final dimensions of the restoration. Instruments with a larger diameter (Step Bur 12S, 1.2 mm) can withstand more milling cycles and have a higher excavation rate than instruments with a smaller diameter (Step Bur 12, 1.0 mm). However, a small diameter is needed to ensure the accurate milling out of smaller and deeper structures. Importantly in this context, the 5-axial milling unit (Arctica; KaVo) incorporates several different instruments, with instrument diameters from 3.0 mm to 0.5 mm. Additionally, with an additional axis, surfaces with an angle close to the insertion axis can be processed more effectively and accurately. The different instruments available within the system make it possible to generate a more accurate relief with deeper fissures and more accurately milled pointed angles.

Provided that quality is not sacrificed, a more rapid milling process for chairside restorations is beneficial for both the dentist and the patient. Changing to smaller instruments and milling with more than 2 different instruments, or milling in several steps, takes time. The milling process with a 5-axial milling unit is slower than that with a 4-axial milling unit, and the 2-step mode requires more time than the 1-step mode. In this context, the faster milling process results in a less accurate restoration, more marginal chipping, reduced longevity, and the possible need to loop in the restoration intraorally.

Future research should investigate the roughness of the edges of restorations milled with different milling strategies (with the scanning electron microscope) and the effect of several mill cycles with the same rotary instrument on milled restorations.

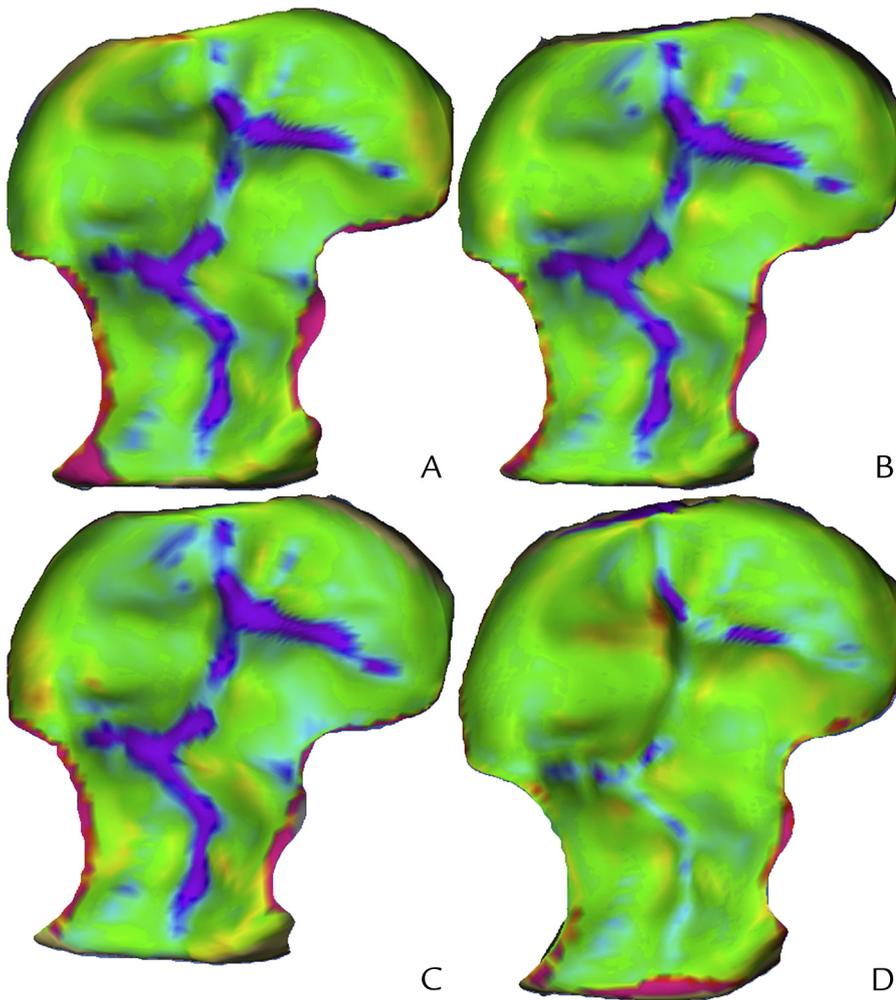

3 Deviations between digital calculated casts and scanned milled casts (trueness). A, Bur 12S, 1-step mode. B, Bur 12, 1-step mode. C, Bur 12, 2-step mode. D, Arctica. Color-coded from $-100\ \mu\text{m}$ (blue) to $+100\ \mu\text{m}$ (red).

The limitations of this research were that only the trueness of the milling process was measured (and not the influence of the discrepancy on the marginal and internal gap) and that the precision of the milling process was not investigated with more milling cycles reusing the same instrument. Scanning with the need to powder the surfaces introduces a certain error level. A scan without powder could reduce possible artifacts.

From the scanning process to the restoration design, the CAD/CAM process has proven accurate under optimal circumstances. The incorporation of an optimized milling process with instruments that have a small diameter may result in similar levels of accuracy and finished restorations with higher aesthetics and better longevity.

CONCLUSIONS

Within the limitations of this *in vitro* study, restorations milled with a 5-axial milling unit exhibited the best quality and the highest trueness values. With regard to the 4-axial milling unit (inLab MC XL; Sirona Dental Systems), the 2-step mode was not significantly better than the 1-step mode. A rotary instrument with a smaller diameter resulted in more accurate milling.

REFERENCES

1. McLaren EA. All-ceramic alternatives to conventional metal-ceramic restorations. *Compend Contin Educ Dent* 1998; 19:307-12.
2. Giordano RA. Dental ceramic restorative systems. *Compend Contin Educ Dent* 1996;17:779-86.

3. Raigrodski AJ. Contemporary all-ceramic fixed partial dentures: a review. *Dent Clin North Am* 2004;48:531-44.
4. Tinschert J, Natt G, Hassenpflug S, Spiekermann H. Status of current CAD/CAM technology in dental medicine. *Int J Comput Dent* 2004;7:25-45.
5. Yeo IS, Yang JH, Lee JB. In vitro marginal fit of three all-ceramic crown systems. *J Prosthet Dent* 2003;90:459-64.
6. Miyazaki T, Hotta Y, Kunii J, Kuriyama S, Tamaki Y. A review of dental CAD/CAM: current status and future perspectives from 20 years of experience. *Dent Mater J* 2009; 28:44-56.
7. Alkumru H, Hullah WR, Marquis PM, Wilson HJ. Factors affecting the fit of porcelain jacket crowns. *Br Dent J* 1988;164:39-43.
8. Balkaya MC, Cinar A, Pamuk S. Influence of firing cycles on the margin distortion of 3 all-ceramic crown systems. *J Prosthet Dent* 2005;93:346-55.
9. Gemalmaz D, Alkumru HN. Marginal fit changes during porcelain firing cycles. *J Prosthet Dent* 1995;73:49-54.
10. Mously HA, Finkelman M, Zandparsa R, Hirayama H. Marginal and internal adaptation of ceramic crown restorations fabricated with CAD/CAM technology and the heat-press technique [published online April 29, 2014]. *J Prosthet Dent*.
11. Bindl A, Mormann WH. Clinical and SEM evaluation of all-ceramic chair-side CAD/CAM-generated partial crowns. *Eur J Oral Sci* 2003;111:163-9.
12. Reiss B, Walther W. Clinical long-term results and 10-year Kaplan-Meier analysis of Cerec restorations. *Int J Comput Dent* 2000;3:9-23.
13. Reiss B. Clinical results of Cerec inlays in a dental practice over a period of 18 years. *Int J Comput Dent* 2006;9:11-22.
14. Luthardt RG, Loos R, Quaas S. Accuracy of intraoral data acquisition in comparison to the conventional impression. *Int J Comput Dent* 2005;8:283-94.
15. Ziegler M. Digital impression taking with reproducibly high precision. *Int J Comput Dent* 2009;12:159-63.
16. Ender A, Mehl A. Full arch scans: conventional versus digital impressions—an *in-vitro* study. *Int J Comput Dent* 2011;14:11-21.
17. Mehl A, Ender A, Mormann W, Attin T. Accuracy testing of a new intraoral 3D camera. *Int J Comput Dent* 2009;12:11-28.
18. Ender A, Mehl A. Influence of scanning strategies on the accuracy of digital intraoral scanning systems. *Int J Comput Dent* 2013;16:11-21.
19. Yuksel E, Zaimoglu A. Influence of marginal fit and cement types on microleakage of all-ceramic crown systems. *Braz Oral Res* 2011;25:261-6.
20. Bindl A, Mormann WH. Marginal and internal fit of all-ceramic CAD/CAM crown-copings on chamfer preparations. *J Oral Rehabil* 2005;32:441-7.
21. Nakamura T, Dei N, Kojima T, Wakabayashi K. Marginal and internal fit of Cerec 3 CAD/CAM all-ceramic crowns. *Int J Prosthodont* 2003;16:244-8.

22. Euán R, Figueras-Álvarez O, Cabratosa-Termes J, Oliver-Parra R. Marginal adaptation of zirconium dioxide copings: influence of the CAD/CAM system and the finish line design [published online January 17, 2014]. *J Prosthet Dent*.
23. Addi S, Hedayati-Khams A, Poya A, Sjogren G. Interface gap size of manually and CAD/CAM-manufactured ceramic inlays/onlays in vitro. *J Dent* 2002;30:53-8.
24. Schaefer O, Watts DC, Sigusch BW, Kuepper H, Guentsch A. Marginal and internal fit of pressed lithium disilicate partial crowns in vitro: a three-dimensional analysis of accuracy and reproducibility. *Dent Mater* 2012;28:320-6.
25. Lee KB, Park CW, Kim KH, Kwon TY. Marginal and internal fit of all-ceramic crowns fabricated with two different CAD/CAM systems. *Dental Mater J* 2008;27:422-6.
26. Boening KW, Wolf BH, Schmidt AE, Kastner K, Walter MH. Clinical fit of Procera AllCeram crowns. *J Prosthet Dent* 2000;84:419-24.
27. Song TJ, Kwon TK, Yang JH, Han JS, Lee JB, Kim SH, et al. Marginal fit of anterior 3-unit fixed partial zirconia restorations using different CAD/CAM systems. *J Adv Prosthodont* 2013;5:219-25.
28. Hamza TA, Ezzat HA, El-Hossary MM, Katamish HA, Shokry TE, Rosenstiel SF. Accuracy of ceramic restorations made with two CAD/CAM systems. *J Prosthet Dent* 2013;109:83-7.
29. Brawek PK, Wolfart S, Endres L, Kirsten A, Reich S. The clinical accuracy of single crowns exclusively fabricated by digital workflow—the comparison of two systems. *Clin Oral Invest* 2013;17:2119-25.
30. Mormann WH, Schug J. Grinding precision and accuracy of fit of CEREC 2 CAD-CIM inlays. *J Am Dent Assoc* 1997;128:47-53.
31. Sjogren G. Marginal and internal fit of four different types of ceramic inlays after luting: an in vitro study. *Acta Odontol Scand* 1995;53:24-8.
32. Keshvad A, Hooshmand T, Asefzadeh F, Khalilnejad F, Alihemmati M, Van Noort R. Marginal gap, internal fit, and fracture load of leucite-reinforced ceramic inlays fabricated by CEREC inLab and hot-pressed techniques. *J Prosthodont* 2011;20:535-40.
33. Felden A, Schmalz G, Hiller KA. Retrospective clinical study and survival analysis on partial ceramic crowns: results up to 7 years. *Clin Oral Investig* 2000;4:199-205.
34. Federlin M, Krifka S, Herpich M, Hiller KA, Schmalz G. Partial ceramic crowns: influence of ceramic thickness, preparation design and luting material on fracture resistance and marginal integrity in vitro. *Oper Dent* 2007;32:251-60.
35. Ilie N, Kunzelmann KH, Hickel R. Evaluation of micro-tensile bond strengths of composite materials in comparison to their polymerization shrinkage. *Dent Mater* 2006;22:593-601.
36. Kramer N, Lohbauer U, Frankenberger R. Adhesive luting of indirect restorations. *Am J Dent* 2000;13:60-76.
37. Lohbauer U, Kramer N, Petschelt A, Frankenberger R. Correlation of in vitro fatigue data and in vivo clinical performance of a glass ceramic material. *Dent Mater* 2008;24:39-44.
38. Groten M, Axmann D, Probster L, Weber H. Determination of the minimum number of marginal gap measurements required for practical in-vitro testing. *J Prosthet Dent* 2000;83:40-9.
39. Moldovan O, Luthardt RG, Corcodel N, Rudolph H. Three-dimensional fit of CAD/CAM-made zirconia copings. *Dent Mater* 2011;27:1273-8.
40. Luthardt RG, Kuhmstedt P, Walter MH. A new method for the computer-aided evaluation of three-dimensional changes in gypsum materials. *Dent Mater* 2003;19:19-24.

Corresponding author:

Gabriel Bosch
Division for Computerized Restorative Dentistry
Center of Dental Medicine
University of Zürich
Plattenstrasse 11
8032 Zürich
SWITZERLAND
E-mail: gabriel.bosch@zsm.uzh.ch

Copyright © 2014 by the Editorial Council for
The Journal of Prosthetic Dentistry.