



**University of
Zurich**^{UZH}

**Zurich Open Repository and
Archive**

University of Zurich
University Library
Strickhofstrasse 39
CH-8057 Zurich
www.zora.uzh.ch

Year: 2012

**Side-by-side comparison of fluoroscopy, 2D and 3D TEE during
percutaneous edge-to-edge mitral valve repair**

Faletra, Francesco F ; Pedrazzini, Giovanni ; Pasotti, Elena ; Moccetti, Tiziano

DOI: <https://doi.org/10.1016/j.jcmg.2012.02.014>

Posted at the Zurich Open Repository and Archive, University of Zurich

ZORA URL: <https://doi.org/10.5167/uzh-75116>

Journal Article

Originally published at:

Faletra, Francesco F; Pedrazzini, Giovanni; Pasotti, Elena; Moccetti, Tiziano (2012). Side-by-side comparison of fluoroscopy, 2D and 3D TEE during percutaneous edge-to-edge mitral valve repair. JACC. Cardiovascular Imaging, 5(6):656-661.

DOI: <https://doi.org/10.1016/j.jcmg.2012.02.014>



Side-by-Side Comparison of Fluoroscopy, 2D and 3D TEE During Percutaneous Edge-to-Edge Mitral Valve Repair

Francesco F. Faletra, MD, Giovanni Pedrazzini, MD, Elena Pasotti, MD, Tiziano Moccetti, MD

CATHETER-BASED MITRAL VALVE CLIP REPAIR IS AN effective procedure in selected patients with mitral regurgitation. Fluoroscopy and 2-dimensional transesophageal echocardiography (2D TEE) are the primary imaging modalities for guidance mitral-clip procedure. Real-time 3-dimensional transesophageal echocardiography (RT 3D TEE) has been recently suggested as a “complementary” imaging modality. However, whether the use of this mitral-clip technique may offer real benefits over standard imaging techniques is not yet established. The following collage, taken from a series of 22 patients, shows side-by-side images obtained by fluoroscopy, 2D TEE, and RT 3D TEE. While essential data during the procedure are provided by the use of the 2 standard imaging techniques and the most crucial step of the procedure (i.e., the capture of mitral leaflets) is always guided by 2D TEE, our collage demonstrates that in almost every step of the mitral-clip procedure, RT 3D TEE may provide new additional useful data. The more precise anatomic information, the fine details of the devices and the precise relationship of catheters/clip with surrounding anatomic structures, may enhance the confidence of imaging interpretation and eventually improve the efficiency of the procedure.

From the Division of Cardiology, Fondazione Cardiocentro Ticino, Lugano, Switzerland. Dr. Faletra has received speaker fees from Philips. All other authors have reported that they have no relationships relevant to the contents of this paper to disclose.

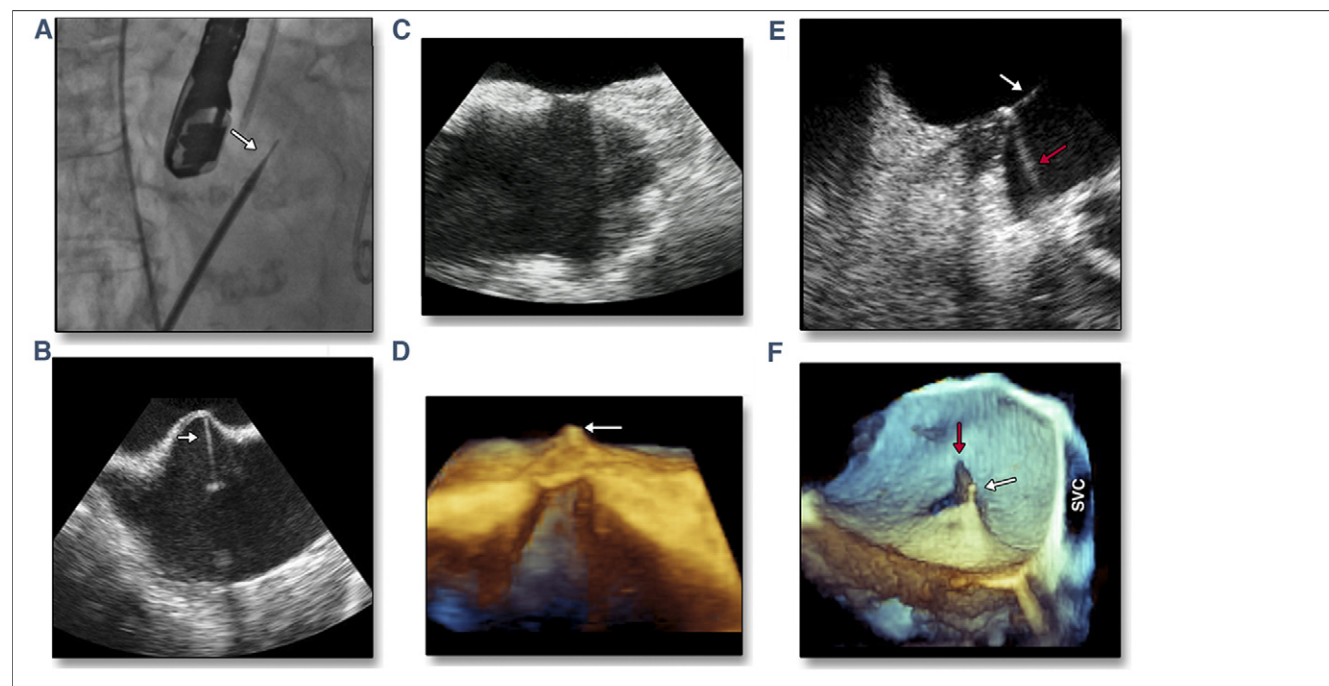


Figure 1. Trans-Septal Crossing

(A) Fluoroscopy of trans-septal crossing. The limited resolution for soft tissue differentiation makes a precise localization of the fossa ovalis difficult (the **arrow** points at the tip of the needle). (B) Two-dimensional transesophageal echocardiography (2D TEE) identification of “tenting”; the **arrow** points at the catheter in right atrium. (C, D) Occasionally, 2D TEE may not show any tenting. Real-time 3-dimensional transesophageal echocardiography (RT 3D TEE) imaging shows that the tenting was deeper in relation to the 2D plane (**arrow**); (E, F) 2D TEE and RT 3D TEE imaging of a deep tenting due to an elastic fossa ovalis. The **white arrows** point to the needle, the **red arrows** point to the reverberations due to the metal. The RT 3D TEE imaging depicts the real anatomic aspect of tenting (“tent-shaped” configuration), which may facilitate understanding of septal morphology and increase the confidence for interventionalists. SVC = superior vena cava.

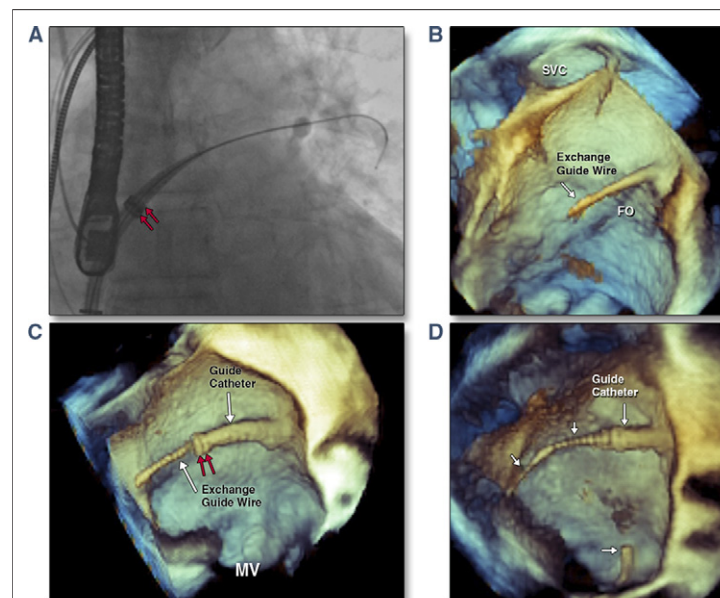
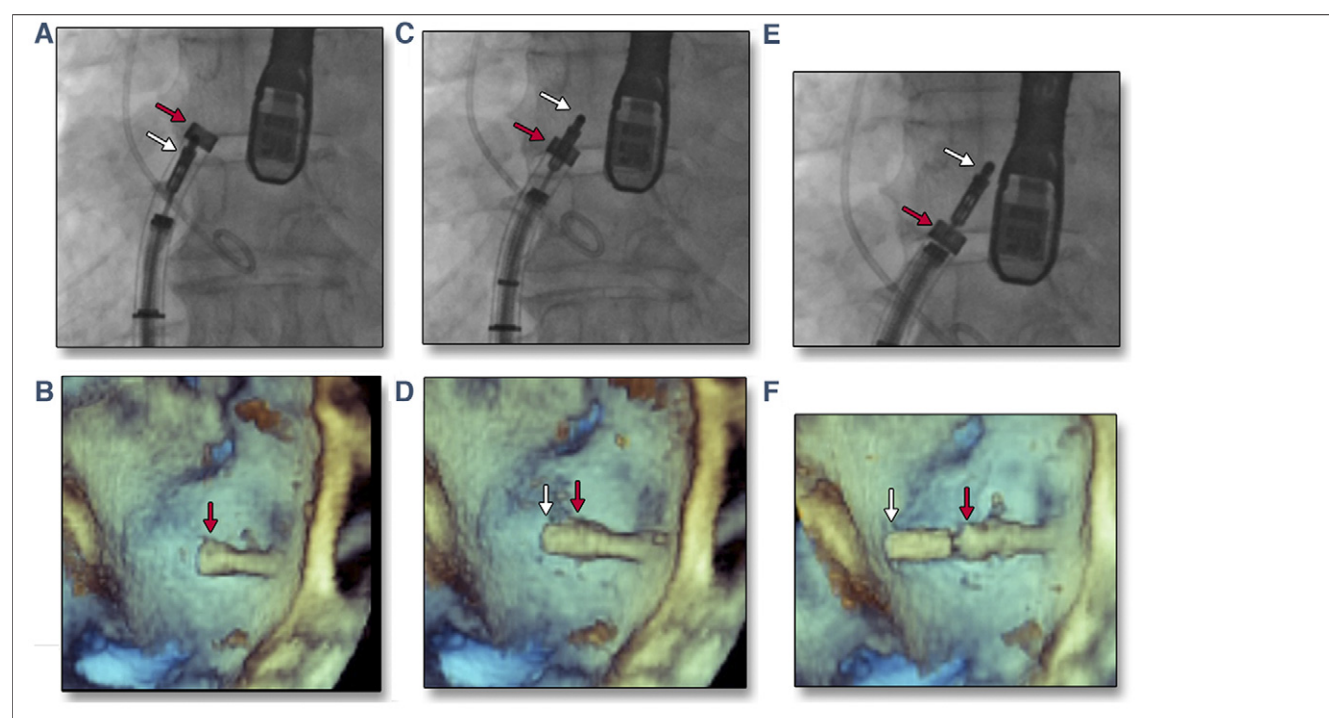
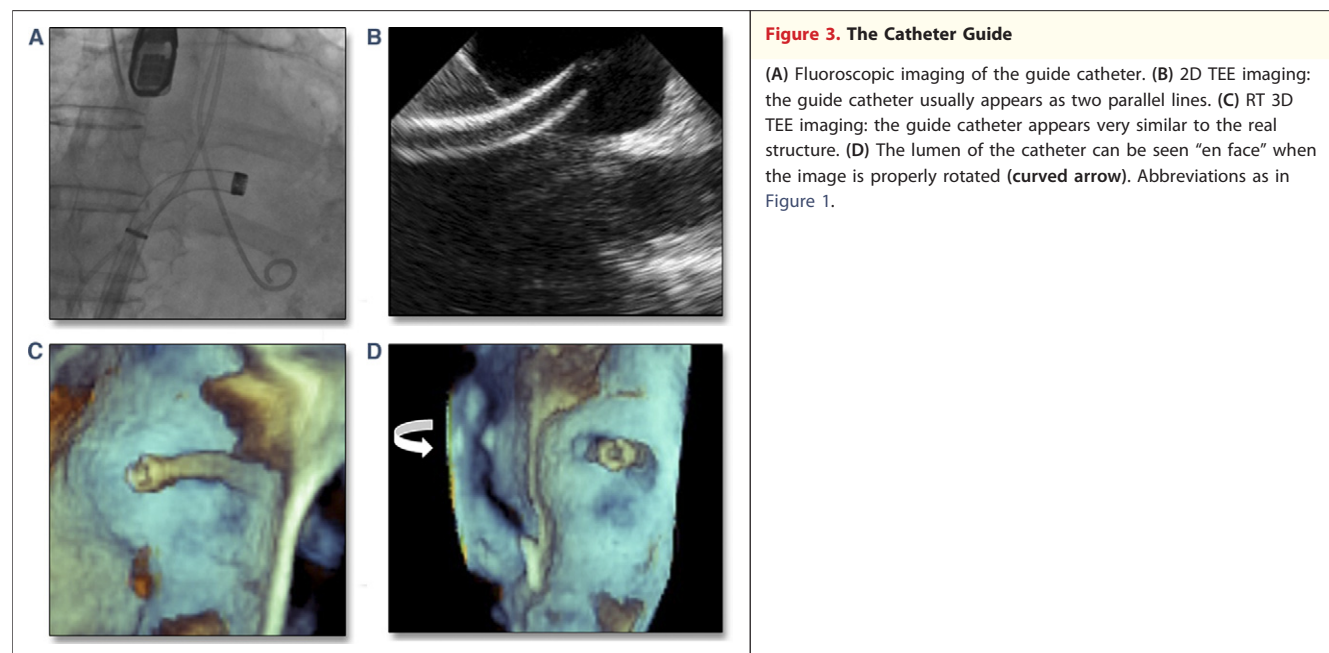


Figure 2. Left Atrial Positioning of Exchange Guide Wire and Guide Catheter

(A) Fluoroscopy imaging of exchange guide wire. The **red arrows** point to the radio-opaque rings used for fluoroscopic identification of the outer catheter guide tip. (B) Exchange guide wire and (C) guide catheter into the left atrium. The radio-opaque rings are easily recognizable (**red arrows**). (D) A long segment of guide wire (**white arrows**) visualized by RT 3D TEE. Online Video 1 shows withdrawal of the guide wire. FO = fossa ovalis; MV = mitral valve; other abbreviations as in Figure 1.



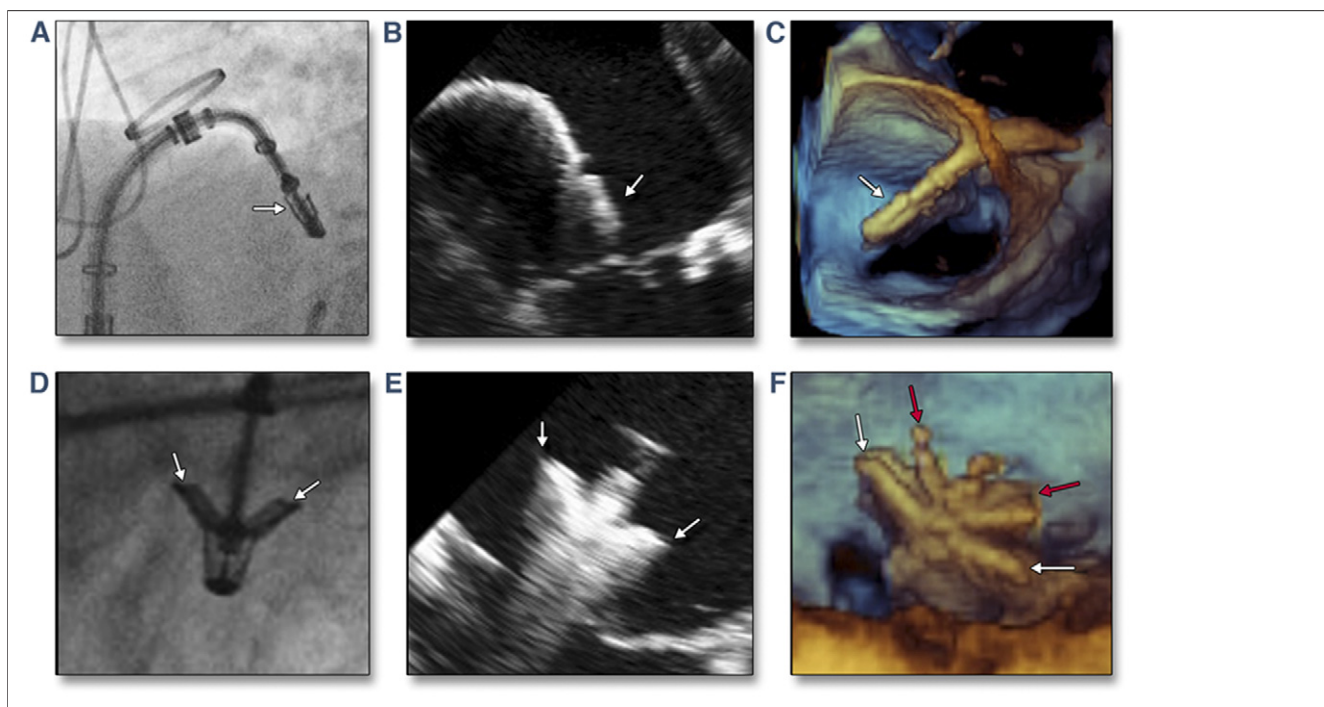


Figure 5. Arms

Arms in closed (A, B, C) and in open (D, E, F) position visualized with (A, D) fluoroscopy, (B, E) 2D TEE, and (C, F) RT 3D TEE imaging. The **arrows** point at the arms in (A, B) closed and (D, E) opened position. RT 3D TEE shows the fine details of arms closed and opened (**white arrows in C and F**) and of grippers (**red arrows in F**). Abbreviations as in Figure 1.

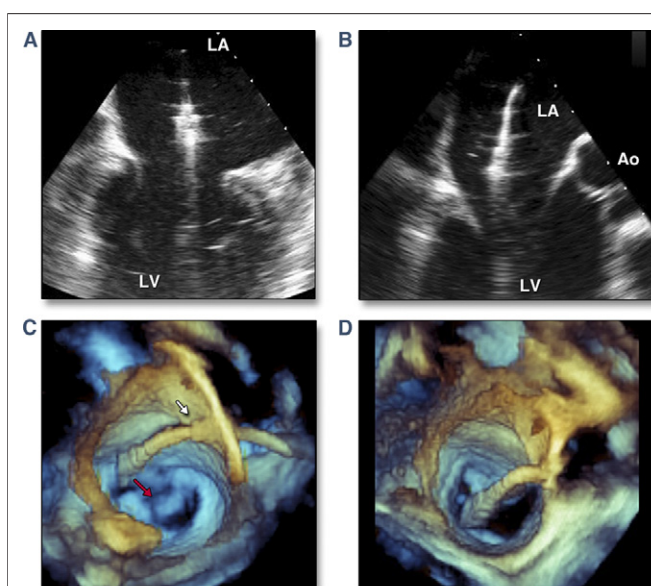


Figure 6. Axial Alignment of the Clip Delivery System

(A, B) 2D TEE guiding axial alignment of the clip. For this maneuver, repetitive cross-checking of two 2-dimensional orthogonal planes (i.e., long axis and intercommissural planes for anterior/posterior and medial/lateral positioning) are required. (C) RT 3D TEE imaging of the mitral-clip system positioned near the lateral commissure. The **white arrow** points to the position of the tip of the guide catheter across the interatrial septum. (D) The mitral-clip system is repositioned (**red arrow in C**) in the center of mitral valve. Ao = aorta; LA = left atrium; LV = left ventricle; other abbreviations as in Figure 1.

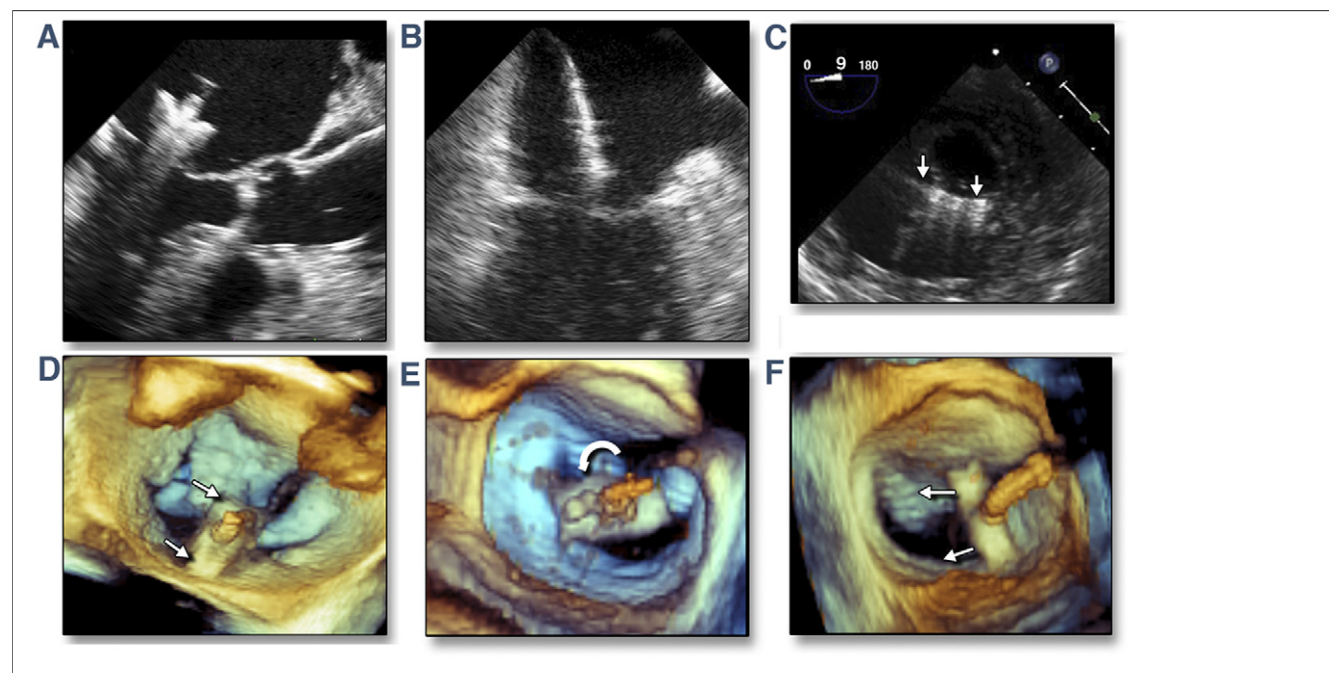


Figure 7. Alignment of the Clip Arms Perpendicular to the Line of Coaptation

(A, B) 2D TEE orthogonal planes must be used for arms alignment. Repetitive cross-checking of these 2 orthogonal planes are usually required. (C) Supplemental transgastric short-axis with the arms perpendicular to the line of coaptation. (D) RT 3D TEE imaging allows an immediate perception of the position of the arms relative to a coaptation line (arrows). (E) Clip arms are not perpendicular to the line of leaflet coaptation. Adjustment requires a counter clock wise rotation (curved arrow); (F) clip arms are oriented perpendicular to the line of leaflet coaptation but positioned near the medial commissure. Repositioning requires a movement toward the center (arrows). Abbreviations as in Figure 1.

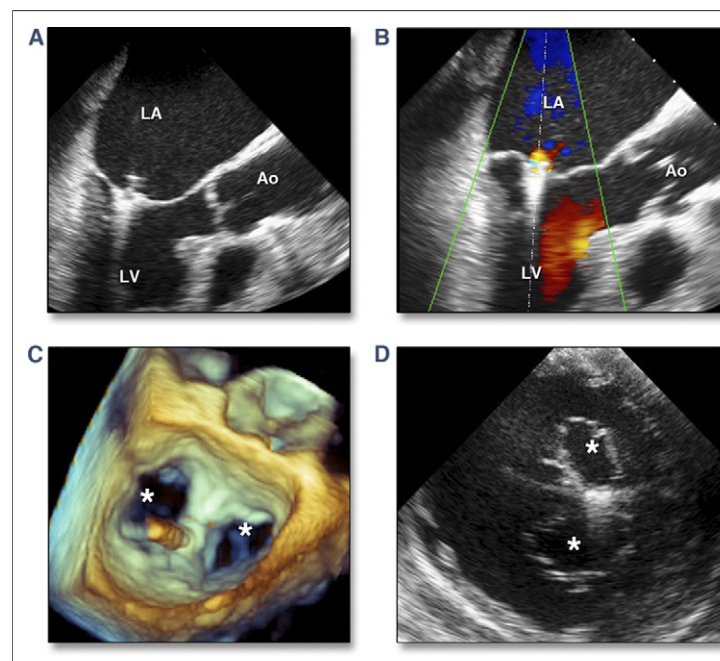


Figure 8. Leaflet Grasping

(A, B) The leaflet capture is guided by 2D TEE. This is because of the higher frame rate and axial resolution in comparison with RT 3D TEE. Moreover, 2D TEE allows the immediate assessment of the residual mitral regurgitation. Although, at present, 2D TEE appears superior to RT 3D TEE, minimizing the pyramidal angle may increase the frame rate up to 30 Hz (Online Video 2). (C) The “double orifice” (*) by RT 3D TEE from the above perspective. (D) The corresponding 2D TEE imaging obtained in short axis transgastric view. Abbreviations as in Figures 1 and 6.

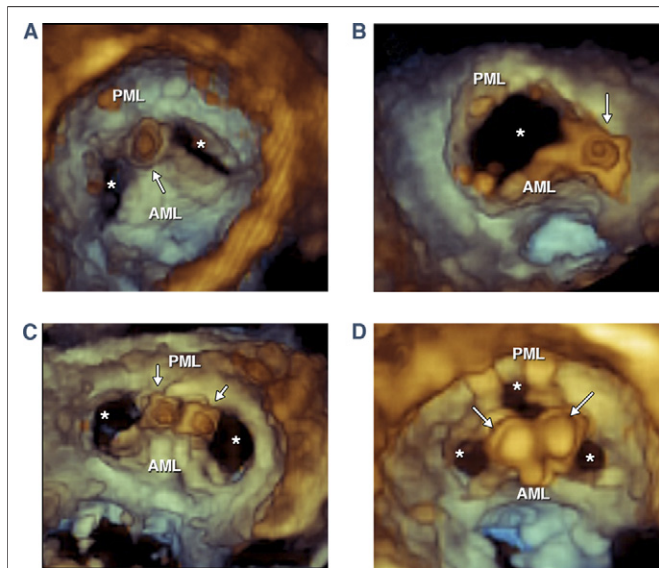


Figure 9. Single, Double, and Triple Orifices

(A, B, C, D) Shows the morphology of the mitral valve after the procedure. The asterisks point to the orifices, arrows point to the clips. (A) One central clip with 2 orifices; (B) 1 medial paracommissural clip with 1 orifice; (C) 2 adjacent central clips with 2 orifices; (D) 2 central separate clips with 3 orifices. AML = anterior mitral leaflet; PML = posterior mitral leaflet.

Address for correspondence: Dr. Francesco F. Faletta, Division of Cardiology, Fondazione Cardiocentro Ticino, Via Tesserete 48, CH-6900 Lugano, Switzerland. *E-mail:* francesco.faletta@cardiocentro.org.

APPENDIX

For supplementary videos and their legends, please see the online version of this article.