



**University of
Zurich**^{UZH}

**Zurich Open Repository and
Archive**

University of Zurich
Main Library
Strickhofstrasse 39
CH-8057 Zurich
www.zora.uzh.ch

Year: 2013

A modified surgical template with dual function

Volpato, Cláudia A M ; Vasconcellos, Diego K ; Garbelotto, Luis Gustavo D ; Manfro, Rafael ; Özcan, Mutlu

DOI: [https://doi.org/10.1016/S0022-3913\(13\)60364-9](https://doi.org/10.1016/S0022-3913(13)60364-9)

Posted at the Zurich Open Repository and Archive, University of Zurich

ZORA URL: <https://doi.org/10.5167/uzh-89935>

Journal Article

Accepted Version

Originally published at:

Volpato, Cláudia A M; Vasconcellos, Diego K; Garbelotto, Luis Gustavo D; Manfro, Rafael; Özcan, Mutlu (2013). A modified surgical template with dual function. *Journal of Prosthetic Dentistry*, 110(3):232-233.

DOI: [https://doi.org/10.1016/S0022-3913\(13\)60364-9](https://doi.org/10.1016/S0022-3913(13)60364-9)

19459

Modified surgical template with dual function: a technical report

Cláudia A. M. Volpato, DDS, PhD,^a Diego K. Vasconcellos, DDS, PhD,^b Luis Gustavo D. Garbelotto, DDS, MSc,^c Rafael Manfro, DDS, MSc^d and Mutlu Özcan, DDS, DMD, PhD^e
School of Dentistry, Federal University of Santa Catarina, Florianópolis, Santa Catarina, Brazil;
North Minas United College/FUNORTE, Florianópolis, Brazil; Center for Dental and Oral
Medicine, University of Zürich, Switzerland

^aAdjunct Professor, Department of Dentistry, Federal University of Santa Catarina.

^bAdjunct Professor, Department of Dentistry, Federal University of Santa Catarina.

^cVisiting Professor, Department of Dentistry, Federal University of Santa Catarina.

^dChairman of Implantodontic Post-graduation, North Minas United College/FUNORTE.

^eProfessor and Head, Dental Materials Unit, Center for Dental and Oral Medicine, University of Zürich.

Corresponding author:

Dr Cláudia Ângela Maziero Volpato

145, Prefeito Antenor Mesquita St.

1104/B

Florianópolis, Santa Catarina

BRAZIL

ZIP: 88015-150

E-mail: claudiamvolpato@hotmail.com, claudiamvolpato@ccs.ufsc.br

19457

A modified surgical template with dual function

ABSTRACT

A modified surgical template to guide the positioning of an implant immediately after tooth extraction and secure the immediate placement of an interim crown is described.

INTRODUCTION

The placement of dental implants immediately after tooth extraction is a popular treatment option, which has been reported to have high success rates.¹⁻⁴ This technique has enabled clinicians to shorten treatment time and possibly provide an interim crown at the time of implant placement. In addition to shortening the surgical procedures, the immediate placement and provisional restoration of a dental implant is purported to better preserve the gingival and bone architecture, resulting in greater comfort and patient acceptance.^{2,4,5}

However, successful immediate implant placement depends on satisfactory patient selection, treatment planning, and atraumatic tooth extraction. This enables adequate primary stability and correct implant location in a site that has preserved the natural architecture of soft and hard tissues. For this purpose, surgical guides have been used to ensure the correct 3-dimensional (3-D) positioning of the implant.^{6,7}

After immediate implant placement, an interim crown can provide patients with rewarding and rapid esthetic results. When an interim crown can be made by copying the

emergence profile of a correctly located but hopeless tooth, it facilitates the maintenance of the postextraction gingival architecture, minimizing possible damage to hard and soft tissues after tooth extraction. However, it is important to ensure that the interim crown does not contact the opposing dentition in maximal intercuspation or during excursive mandibular movements.^{9,10} In this article, a surgical guide that can correctly guide implant placement as well as the correct placement of an interim crown is described.

TECHNIQUE

1. Make stone casts of the maxillary and mandibular arches. Remove the tooth to be extracted from the cast and make an interim crown from heat-polymerized polymethyl methacrylate (New Outline; AnaxDent GmbH, Stuttgart, Germany).
2. Fabricate a clear plastic template in 2-mm thick, clear, soft ethylene vinyl acetate (Pro-Form; Dental Resources Inc, Delano, Minn) on the cast with the interim crown (Fig 1). Cut the acetate sheet following the buccal gingival contour. Keep the guide open on the palatal side around the crown to enable surgical access.
3. Extract the tooth by using an atraumatic technique (Benex - Control Root Extraction System; Meisinger, Düsseldorf, Germany) (Fig 2). After extraction, place the surgical guide to direct the positioning of the implant. The open palatal side of the guide will enable surgical access, and the buccal outline will assist proper 3-D implant positioning (Fig 3).
4. After implant placement (Flash 3.5 × 13mm; Conexão, São Paulo, Brazil), attach the corresponding provisional abutment to the implant to fabricate the interim crown. Evaluate whether adequate occlusal space exists for the crown and the abutment. If not, unscrew the component (Abutment 166211; Conexão) and shorten it outside the mouth. Condense a small

piece (approximately 5 mm) of polytetrafluoroethylene tape (PTFE; New Age Industries Inc, Willow Grove, Pa) in the screw hole to prevent acrylic resin from sticking to the head of the screw.

5. Place the interim crown in the previously fabricated surgical guide. With a brush-bead technique (Figs 4 and 5), fill the crown with colored autopolymerizing acrylic resin (Duralay; Reliance, Dental MFG Co., Worth, Ill). Position the assembly intraorally. Keep the surgical guide in place until the acrylic resin polymerizes. Then, remove the acetate guide, access the abutment screw with a round diamond rotary instrument (Komet, Rock Hill, SC), remove the crown and create the appropriate emergence profile with acrylic resin (Duralay; Reliance, Dental MFG. Co.).

6. Replace the crown to verify whether the gingival architecture is preserved. In the event of excessive pressure, reduce the emergence profile. Finish and polish (Acrylic Polisher; Komet) the crown, especially in the transmucosal area (Fig 6).

SUMMARY

This technique allows the clinician to use a single surgical guide to correctly place an implant in a postextraction site and position an interim crown that restores optimal esthetics.

REFERENCES

1. Groisman M, Frossard WM, Ferreira HM, De Menezes Filho LM, Touati B. Single-tooth implants in the maxillary incisor region with immediate provisionalization: 2-year prospective study. *Pract Proced Aesthet Dent* 2003;15:115-22.
2. Kan JY, Rungcharassaeng K, Lozada J. Immediate placement and provisionalization of maxillary anterior single implants: 1-year prospective study. *Int J Oral Maxillofac Implants* 2003;18:31-9.
3. Covani V, Crespi R, Cornelini R, Barone A. Immediate implants supporting single-crown restoration: a 4 year prospective study. *J Periodontol* 2004;75:982-8.
4. Cosyn J, Eghbali A, De Bruyn H, Collys K, Cleymaet R, De Rouck T. Immediate single-tooth implants in the anterior maxilla: 3-year results of a case series on hard and soft tissue response and aesthetics. *J Periodontol* 2011;38:746-53.
5. Den Hartog L, Raghoobar GM, Stellingsma K, Meijer HJ. Immediate loading and customized restoration of a single implant in the maxillary esthetic zone: a clinical report. *J Prosthet Dent* 2009;102:211-5.
6. Schneider D, Marquardt P, Zwahlen M, Jung RE. A systematic review on the accuracy and the clinical outcome of computer-guided template-based implant dentistry. *Clin Oral Implants Res* 2009;4:73-86.
7. Oh TJ, Shotwell J, Billy E, Byun HY, Wang HL. Flapless implant surgery in the esthetic region: advantages and precautions. *Int J Periodontics Restorative Dent* 2007;27:27-33.
8. Al-Harbi SH. Minimizing trauma during tooth removal: a systematic sectioning approach. *Eur J Esthet Dent* 2010;5:274-87.

9. Valentini P, Abensur D, Albertini JF, Rocchesani M. Immediate provisionalization of single extraction-site implants in the esthetic zone: a clinical evaluation. *Int J Periodontics Restorative Dent* 2010;30:41-51.
10. Balshi TJ, Wolfinger GJ, Wulc D, Balshi SF. A prospective analysis of immediate provisionalization of single implants. *J Prosthodont* 2011;20:10-5.

LEGENDS

Fig. 1. Acetate sheet on cast with interim crown.



Fig. 2. Atraumatic tooth removal.



Fig. 3. Surgery assisted with guide.



Fig. 4. With brush-bead technique interim crown is filled with autopolymerizing acrylic resin.

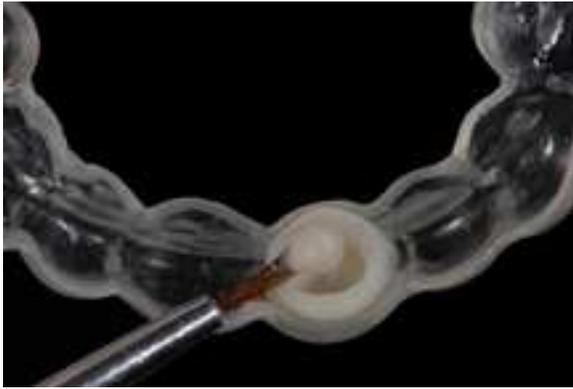


Fig. 5. Guide and interim crown intraorally.



Fig. 6. Result after placement of interim crown.

